

*A summary of new products and services for materials research . . .*

**Rare Earth Materials:** Megon A/S of Norway is now offering its rare earth products in the United States through Crystalox, Ltd. Megon supplies materials for the growth of laser and microwave crystals and for producing ceramic parts for microwave systems. Employing a patented process, Megon can meet individual customer's specifications from 4N (99.99%) to 6N (99.9999%) purity. Other technical services include custom process development for purification of rare earth mineral resources and concentrates, or material recovery. Crystalox, Ltd., Suite 101, 100 Brush Creek Rd., Santa Rosa, CA 95404-2709; (707) 539-2508.

**LaGaO<sub>3</sub> Crystal Substrates:** Substrates for superconducting thin films are 1.6 cm diameter with (100) orientation. Larger sizes and different materials will be available later. The company's product line also includes crystals and melt cast samples of: Y<sub>1</sub>Ba<sub>2</sub>Cu<sub>3</sub>O<sub>x</sub>, Bi<sub>2</sub>Ca<sub>1</sub>Sr<sub>2</sub>Cu<sub>2</sub>O<sub>x</sub>, CuO, Cu<sub>2</sub>O, and long, thin superconducting crystals of proprietary composition. SuperconiX Inc., 261 East 5th St., Saint Paul, MN 55101; (612) 222-0046.

**SQUID-Based Magnetometer:** SQUID-based variable temperature magnetometer measures the dc magnetic susceptibility of small samples in a wide range of applied magnetic fields and can also be configured for other measurements, including resistivity and Hall Effect. The MPMS<sub>2</sub> magnetometer uses a low T<sub>c</sub> SQUID amplifier to measure magnetic moment over a dynamic range of 10<sup>7</sup> to 1.25 emu. With a differential sensitivity of 10<sup>-6</sup> emu, the MPMS<sub>2</sub> can characterize samples in both thin film and bulk forms up to 8 mm diameter and 7 mm long. During measurement, materials can be exposed to applied magnetic fields, up to +/-10,000 Oe, and temperatures from 5 to 350 K. The MPMS<sub>2</sub> can be easily programmed to run standard tests completely unattended, and can be used for characterizing high T<sub>c</sub> superconductors, spin-glass materials, and hard or soft magnetic materials. The competitively priced system can be configured with several options which allow it to measure permanent or ferromagnetic materials (to +/-300 emu), quickly collect magnetic hysteresis data, or measure resistivity and magnetoresistance. Quantum Design, 11578 Sorrento Valley Rd., Suite 30, San Diego, CA 92121; (619) 481-4400.

**Ultrathin Flexible Zirconia Ceramic:** Commercially available ultrathin zirconium oxide material formed into sheets 0.05 to 0.10 mm thick features high strength and toughness, high resilience and flexibility, and ease of fabrication. The 0.05 mm sheets can be bent around a radius of 5 mm and will return to their original shape when released. CERAFLEX sheets can be formed into rings, spring shapes, corrugated sheets, short tubes, and other curved configurations. Holes can be machined as close as 0.2 mm apart with diameters from 0.01 to 4 mm. The materials is easily laser machined with automatic processing equipment. Applications include superconductor substrates, fuel cells, wafers and substrates for hybrid circuits, and parts for precision instruments and electronics. Marketech International, 414 South Craig St., Suite 300, Pittsburgh, PA 15213; (412) 3103.

**High Pressure Oxygen System for Superconductors:** Self-contained bench-top unit has been developed for the synthesis and high temperature treatment of high temperature superconductors in high pressure oxygen. Along with the ability to change samples rapidly and to change pressure during system operation, features include programmable multi-step intelligent temperature control, and accessories capable of extending the range of measurements and applications. No special installation, utilities, or support services are necessary because the system is a standard oxygen cylinder and is powered by a standard electrical wall outlet. Maximum pressure available is 200 atm of pure oxygen at temperatures up to 1,000°C. Models compatible with U.S., European, and Japanese electric supplies are available. Accessories under development include a quench option for rapid cooling, high pressure oxygen flow-through system for continuous purge, and sample resistivity vs. temperature and oxygen pressure. Morris Research, Inc., 1862 Euclid Ave. #240, Berkeley, CA 94709; (415) 486-4142.

**Solid Arsenic Charges for Ion Implantation:** Solid arsenic charges designed to increase the efficiency of ion implantation manufacturing and operations are available in sizes to fit all major ion implantation units. A one-piece design simplifies handling of toxic material and eliminates difficult, time-consuming vaporizer loading. The high-density charges have a minimum exposed surface area, virtually eliminating surface oxidation and contamination prob-

lems. All charges are individually vacuum packaged in laminated poly/foil bags for easy storage and loading. Johnson Matthey, AESAR Group, Eagles Landing, P.O. Box 1087, Seabrook, NH 03874; (800) 343-1990, in New Hampshire (603) 474-5511.

**Ceramic Abstracts Online:** CERAB, the American Ceramic Society's online version of Ceramic Abstracts, has been added to STN International, a scientific and technical information network operated in North America by Chemical Abstracts Service. CERAB contains more than 100,000 bibliographic records from books, conference proceedings, journals, patents, reports, and trade literature. Coverage spans worldwide scientific, commercial, and engineering literature on ceramics, glass processing, magnetics, materials, superconductors and many other ceramics-related topics. Materials can be searched by chemical formulas, and the 24-hour service covers Ceramic Abstracts since 1976. Reference Services Department, American Ceramic Society, 757 Brooksedge Plaza Drive, Westerville, OH 43081-6136; (614) 890-4700.

**Directory of Ceramics Products and Services:** *Ceramic Source '90*, an annual reference guide for the ceramics field, provides comprehensive alphabetical listings of products, services, product tradenames, and technical data, and an international company directory of over 3,500 listings. The product directory covers raw materials/ceramic powders, production equipment and supplies, coatings and finishings, testing and evaluation, and ceramic components and devices. The technical data section provides tables on ceramic properties published over the last several years. \$28 per copy. American Ceramic Society, 757 Brooksedge Plaza Drive, Westerville, OH 43081-6136; (614) 890-4700.

**Optics/Optoelectronic Network:** New electronic networking service is designed to link optical and optoelectronic scientists and engineers worldwide with a full menu of services and databases offered by SPIE, The International Society for Optical Engineering. OPTO-LINK includes an electronic mail service, the Society's Working Group bulletin boards, membership directory, publications listings, and conference and technical exhibition schedules as well as resume files and job postings. A major expansion is planned in 1990. For access, dial (206) 733-2998 via modem from any computer. SPIE, P.O. Box 10, Bellingham, WA 98227-0010; (206) 676-3290. □



# TEM Specimen Preparation

Urgent materials problem?  
Answer required tomorrow?  
Gatan's range of specimen preparation instruments is designed to produce high quality TEM specimens quickly, efficiently and reliably. We are confident that our equipment will satisfy your needs because we use it ourselves every day in our own commercial specimen preparation service.

## Disc Fabrication

Preparation of good TEM specimens from non-biological materials starts with the production of high quality discs. For this step Gatan offers an ultrasonic disc cutter (hard materials) and a disc punch (soft materials).

## Pre-thinning

To obtain large electron transparent areas the starting discs should be mechanically pre-thinned. For this delicate work, Gatan provides a high precision disc grinder and a dimple grinder. Our dimple grinder is so sensitive that electron transparent specimens can be prepared directly.



## DuoMill™

A high flux, broad beam, twin station ion mill with our famous Whisperlok (US Patent 4,272,682) specimen exchange system and long life Octoguns™ for high speed final thinning of specimen discs. Iodine/Argon Octoguns™ are now available for reactive ion milling of semiconductors.

## PIMS™

A precision ion milling system (shown here) to thin selected areas at interfaces in semiconductors and other materials. It produces an image to show you where and how much ion milling has occurred. The <math><2\ \mu\text{m}</math> diameter focused ion beam can also be computer controlled to cut any surface pattern you desire.

Please visit Booth No. 605 at the MRS Show in San Francisco, April 17-19, 1990.

For further information on TEM specimen preparation equipment, TV-systems, EELS, PEELS™ and specimen holders, please call our main sales and service offices at:

**Gatan Inc.**, Pleasanton, CA  
Tel. (415) 463-0200

**Gatan, Inc.**, Warrendale, PA  
Tel. (412) 776-5260

**Gatan GmbH**, München, FRG  
Tel. (089) 352374

**Gatan Ltd.**, High Wycombe, UK  
Tel. (0494) 26211


# The **ULTIMATE** Vacuum

VG Special Systems Ltd., the world-leader in ultra high vacuum engineering, builds UHV systems for thin films research, surface analysis and specialized production environments.


If your application demands guaranteed pressures in the mid  $10^{-11}$  torr range and bakeout capability to 450°C call us.



Production UHV for heat treatment

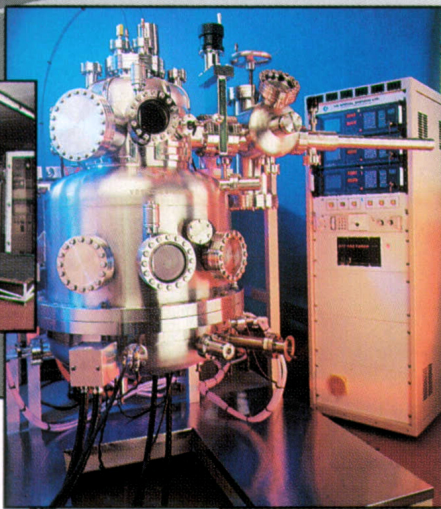


Production UHV for MBE thin films

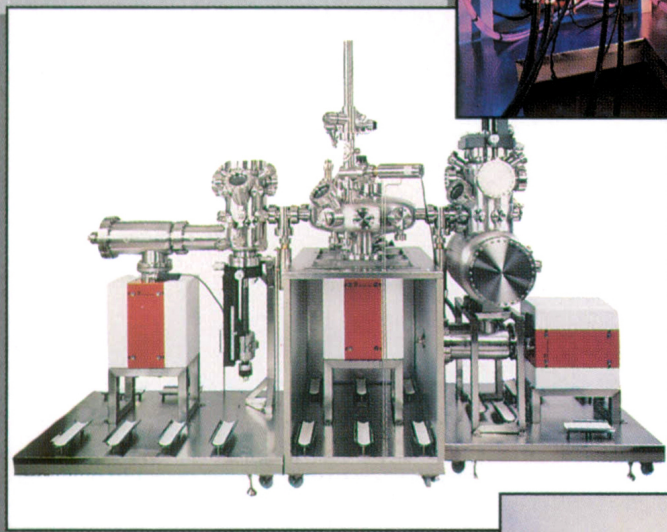


Analytical UHV for multiple surface science techniques



Production UHV for electronics components

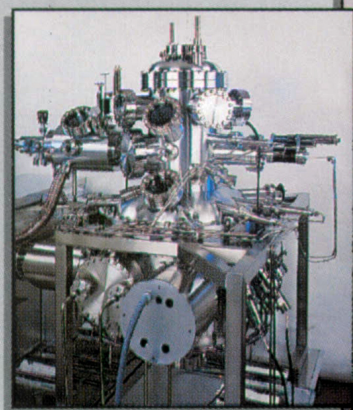

Research UHV for thin film high  $T_c$  superconductors

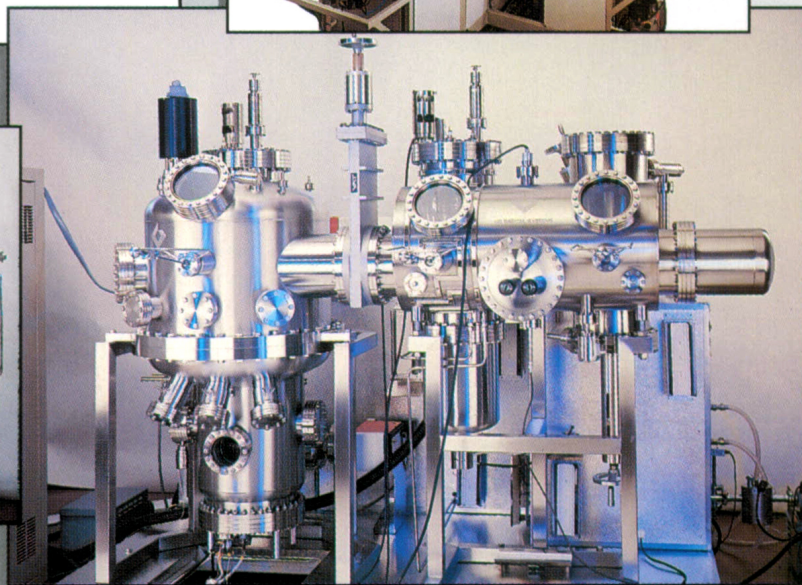

Research UHV for evaporation onto heated 6" wafers



**Kurt J. Lesker**  
Company


1515 Worthington Ave. Clairton PA 15025

"... vacuum science is our business"  
**800-245-1656 USA**

800-242-0599 (PA) 412-233-4200  
800-848-5386 (CA) Fax: 412-233-4275

Please visit Booth No. 414 at the MRS Show in San Francisco, April 17-19, 1990.

*Do you need to position samples to 0.001 mm?  
Cycle samples between 20 and 1300 Kelvin?  
At a base pressure of  $5 \times 10^{-11}$  torr?*

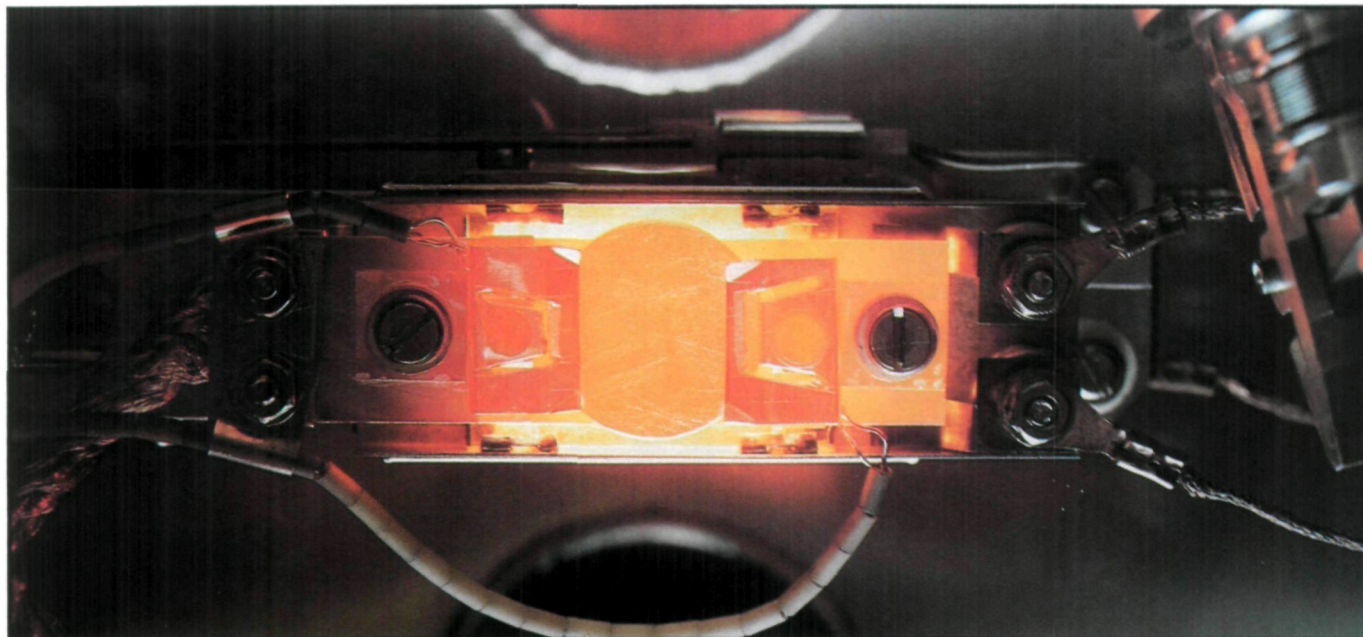


*An LN<sub>2</sub> cooled HPT manipulator operating at high temperature*

*... Then only our manipulators will work!*

Vacuum Generators high precision sample manipulators are designed and constructed with the skill, care, and quality only the world's premier vacuum manufacturer can provide.

- ♦ Short and long travel (up to 600 mm)
- ♦ Wide bore travel (up to 25 mm)
- ♦ Resistance or e-beam heating
- ♦ Resolution 0.0005 mm (on x and y axis)
- ♦ Bakeable to 230°C without dismantling
- ♦ 5-axis and tilt (HPT and HPLT models)
- ♦ 6-axis (CPTO models – rotational resolution 0.005°)
- ♦ Motor drives for all axes
- ♦ LN<sub>2</sub> or LHe cooled versions

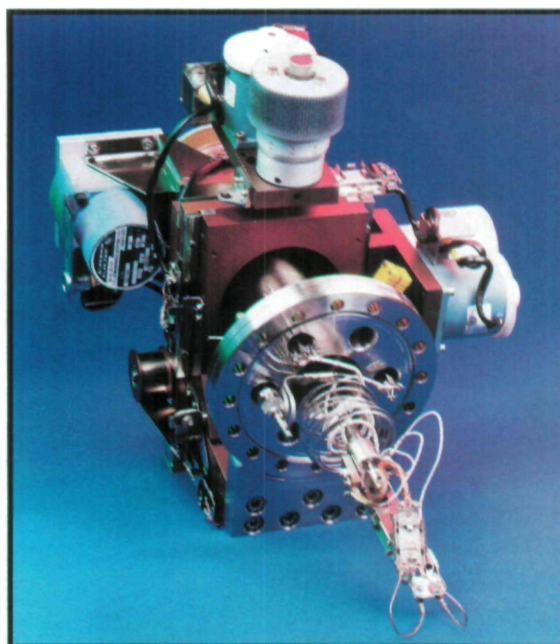

*Fully motorized 5-axis HPT high precision manipulator*

Please visit Booth No. 414 at the MRS Show in San Francisco, April 17-19, 1990.

Vacuum Generators represented in the USA by:



**Kurt J. Lesker**  
Company


1515 Worthington Avenue Clairton, PA 15025 412-233-4200

"... vacuum science is our business"

800-245-1656 USA  
800-848-5386 (CA)

800-544-3940 CAN  
800-242-0599 (PA)

# Measuring Magnetic Susceptibility Has Never Been This Easy— Or This Accurate.

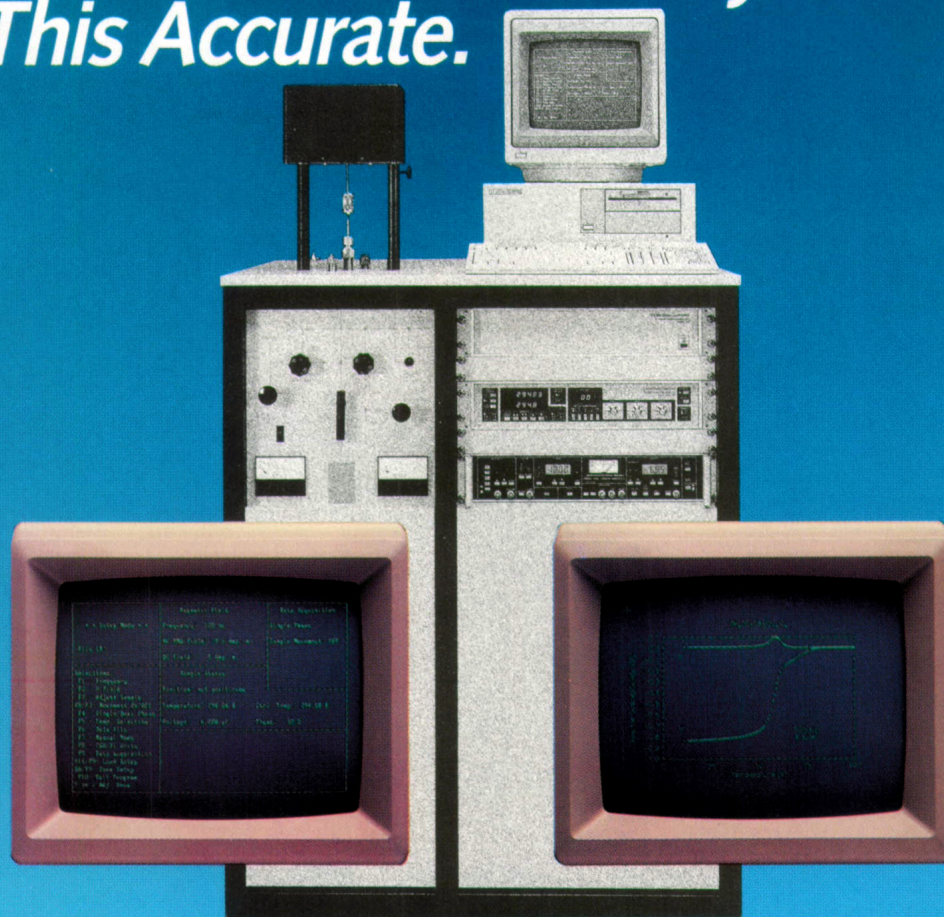

**The Lake Shore Model 7000 AC Susceptometer.**  
High Sensitivity—better than  $10^{-7}$  emu (in magnetic moment)—for measurements at low ac magnetic fields.

## Features:

- Accurate measurements, in absolute values, of both the real ( $X'$ ) and imaginary ( $X''$ ) components of the magnetic susceptibility of ferromagnetic, paramagnetic and diamagnetic materials, including bulk and thin film superconductors.
- Operational flexibility—absolute determinations of ac susceptibility as a function of temperature (<4.2K to 330K), frequency (5 to 1000Hz) and field (0.005 to 10 Oe RMS); resistivity measurements;

customized sample profiling—including temperature sweep or fixed point measurement techniques.

- Computer-automated data acquisition/analysis software.
- Rapid sample throughput.

For accurate characterization of conventional or high temperature superconductors and other magnetic property materials—plus high value, low cost and prompt delivery—discover the Model 7000 AC Susceptometer from Lake Shore.

Please visit Booth No. 505 at the MRS Show in San Francisco, April 17-19, 1990.

Look to  
Lake Shore.  
High-  
performance  
in low  
temperature  
technology

## LakeShore



CRYOTRONICS, INC.

64 East Walnut Street, Westerville, Ohio 43081 USA (614) 891-2243 Telex: 24-5415 Cryotron WTVL Fax: (614) 891-1392

Cryophysics - Oxford, UK; Darmstadt, W. Germany; Jouy-En-Josas, France; Geneva, Switzerland; J.C. Electronics Corp. Tokyo, Japan;  
Goodwill Cryogenics Enterprises - New Bombay, India; SST Banin Enterprise Co., Ltd. - Taipei, Taiwan R.O.C.;  
Woo Sin Enterprise, Inc. - Seoul, Korea; US-China Industrial Exchange - Beijing, P.R.C.

© 1989 Lake Shore Cryotronics